



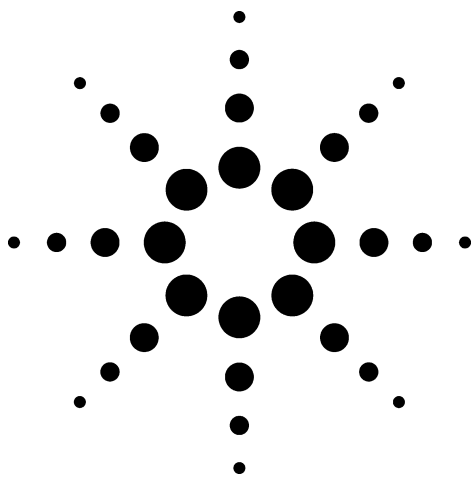

Calmodulin Tryptic Digest Flow Rate Effects

Application
Biochemical
Robert Ricker

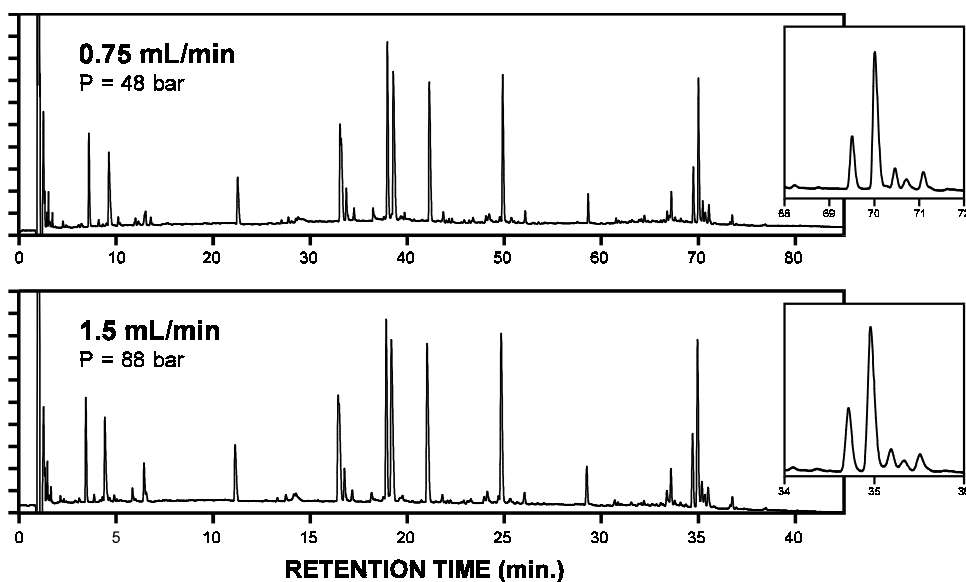

Highlights

- ZORBAX SB-C18 (3.5 μ m) provides high resolution at reasonable back pressures.
- SB-C18 bonded phase is very stable at temperatures up to 90°C and at low pH (e.g. 2.0).

Conditions:
ZORBAX SB-C18 (3.5 μ m), 4.6 x 150 mm (Agilent P/N: 863953-902)
85°C, 45 μ l injection (75 μ g) UV 210 nm
Gradient: 5-60% B in 45 min, (1.5 mL/min) or 90 min (0.75 mL/min)
A-0.1% TFA in water, B-0.095% TFA in 80% ACN/water



Agilent Technologies

*Robert Ricker is an application chemist
based at Agilent Technologies, Wilmington,
Delaware.*

For more information on our products and
services, visit our website at:
www.agilent.com/chem

Copyright© 2002 Agilent Technologies, Inc.
All Rights Reserved. Reproduction,
adaptation or translation without prior
written permission is prohibited, except as
allowed under the copyright laws.

Agilent shall not be liable for errors
contained herein or for incidental or
consequential damages in connection with
the furnishing, performance, or use of this
material.

Information, descriptions, and specifications
in this publication are subject to change
without notice.

Printed in the USA
April 25, 2002
5988-6308EN



Agilent Technologies